HUMAN CARE

Your Life. Your Way.

HeliQ, HeliQ power traverse
55100, 55110, 55120, 55130, 55140, 55150

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Content

ENGLISH

| | |
|---------------------------------------|----|
| User Manual | 21 |
| Compliances and Standards | 21 |
| Safety instructions | 22 |
| Specifications | 23 |
| Product description | 24 |
| Intended use | 24 |
| Product label | 25 |
| Symbols on the unit and product label | 25 |
| Compatibility | 26 |
| Functions | 26 |
| Suspending the lift from the rail | 27 |

| | |
|--|----|
| Hand control | 28 |
| Secondary control panel | 28 |
| Indicator | 29 |
| Charger | 29 |
| Detaching the lift from the rail system | 30 |
| Using the lift | 30 |
| Data communication | 30 |
| Hanger Bars | 31 |
| Lifting height | 31 |
| Accessories | 31 |
| Troubleshooting | 32 |
| Care instructions for the lift and accessories | 33 |
| Guidance and manufacturer's declaration | 34 |

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User Manual

Important!

You must read the user manual for your device prior to usage. Keep this booklet and information handy for future use.

“WARNING!” This symbol is used when particular attention is required

Lifting and moving a person always entails some risk. You should therefore read the user manuals for the lift and the lifting accessories carefully. Always make sure that the accessories you are using are intended for use with the lift. As a caregiver, you are responsible for the client's safety and you must know whether or not the client will be able to cope with the lifting process.

Human Care's lifts are designed only to raise and lower clients in a vertical direction, and never diagonally or through angles. The lift must always be positioned directly above the client when lifting. The client must always be positioned above the surface onto which they are to be lowered.

LIFTING OR LOWERING A Client DIAGONALLY OR IN AN ANGLED MOVEMENT MAY RESULT IN INJURY TO THE Client OR THE CAREGIVER.

If the product is used in any way that is not recommended by Human Care, then Human Care will not accept responsibility for any injuries or accidents that may occur.

If no other agreement or approval exists, the owner of the product is responsible for any combinations using accessories or products from suppliers other than Human Care.

Always contact your supplier if you have any questions or if anything is unclear.

Anyone using the equipment must have been given adequate training and be certified by Human Care in the use of the lift and its accessories

Clients must not be left unattended and should have a trained individual with them when performing a lift.

Human Care's products are continuously being developed and updated and we reserve the right to make design changes without prior notice.

Notice to user/patient in case of serious incident

Any serious incident that has occurred in relation to the product, should be reported to the local contact, who reports to the manufacturer, and the competent authority of the country in which the user/patient is established.

Compliances and Standards

Human Care is an ISO 13485:2016 certified Swedish medical device company. The Quality Management System is in compliance with US 21 CFR part 820.

The product is CE marked in accordance with EU

Directive 93/42/EEC (MDD) and/or EU Regulation 2017/745 (MDR), as class I medical device.

The CE mark is on the product.

This product has been tested in accordance with ISO 10535.

Safety instructions

Before using for the first time, make sure that:

- You have read and understood the user manual for the lift and all of its accessories
- The lift has been installed in accordance with the installation instructions
- The lift accessories are intended for use with the lift and are correctly fitted
- Before use: charge battery upon arrival. To ensure battery capacity is maintained, the new ceiling lift must be charged for 8 hours immediately upon arrival before it is put in use or stored. This is because the ceiling lift may have been stored for some time and during that period not been charged regularly.
- Remove the lift and charger from the packaging.
- Remove the seal from the lift belt after first reading the user manual.
- Plug in the charger at an easily accessible wall socket.
- Plug the hand control into the charger (see page 8). If the indicator light on the charger is glowing yellow, charge until the indicator light turns green.

Important! The lift is designed exclusively for lifting people, and only one at a time.

Before use, always ensure that:

- Personnel who use the equipment have been given adequate training and been certified by Human Care in the use of the lift and its accessories
- The lift belt is not twisted or worn and that it can move freely in and out of the lift
- Lift accessories, e.g. sling, distance strap and rail roller, are showing no signs of wear or other damage

- The lift is correctly suspended in the distance strap, rail roller and rail
- The right types, sizes, materials and designs of the accessories have been selected with respect to the client's safety and needs
- The lift accessories have been correctly installed. The sling's strap loops and suspension must be checked once the loops are fully extended when lifting but before the client is lifted from the underlying surface.
- It is important to always check that the loops are at the bottom of all hooks (see picture below), securely placed below the locking pin.

Right

Wrong

Human Care provides suitable training in safe operation of the lift and other accessories.

Further information about slings, rail systems and other accessories can be found at www.humancaregroup.com.

Important! To ensure that the lift belt will not be damaged, it is important that the lift bar is balanced and that the sling is suspended correctly. Never hold or pull the lift belt when moving a client. It is also important that the lift belt is tensioned as it is wound in.


Warning! Modification of the product is not permitted without the manufacturer's consent!

Specifications

| Lifting capacity, max.: | |
|--------------------------------|---|
| 55100 | 150 kg/330 lbs |
| 55110 | 220 kg/484 lbs |
| 55120 | 300 kg/661 lbs |
| 55130 | 150 kg/330 lbs power traverse |
| 55140 | 220 kg/484 lbs power traverse |
| 55150 | 300 kg/661 lbs power traverse |
| Dimensions: | |
| (WxHxL) | 235x140x235mm / 9.2x5.5x9.2" |
| Weight of lift: | 7 kg/15.4 lbs |
| Batteries: | 2x12V, 2.3Ah |
| Emergency lowering: | Electronic & Mechanical |
| Class: | Internally powered equipment |
| Type: | B |
| IP rating lift: | IPX4 |
| IP rating hand control: | IPX7 |
| Noise level: | <65dB (loaded and unloaded) |
| Lifting speed: | 1.8 m/min / 3.0 m/min (unloaded only) |
| Lifting height: | 2120 mm/ 83,5" (150/220 kg), 10680/66,1" (300 kg) |
| Start/stop: | soft |
| Indicator: | |
| Standby | flashing GREEN |
| Active | constant GREEN |
| Error | constant RED |
| Service | flashing YELLOW/GREEN, YELLOW/RED |
| Charging | constant YELLOW |
| Charger | Input 100-240V AC 50-60 Hz, 1.6A; Output 27.8V, 0.8A DC |
| IP rating charger: | IPX0 |
| ETL: | Yes |

This product is designed for repeated use.

Product description

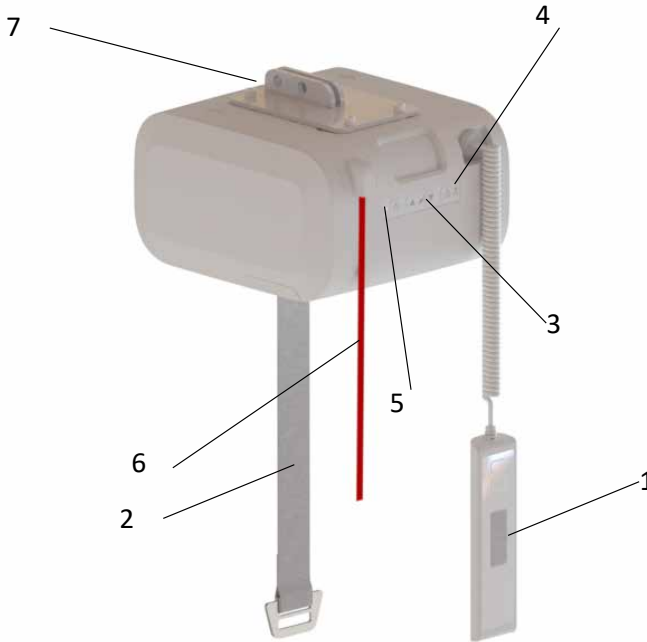1. Hand control
2. Lift belt
3. Secondary control panel
4. Emergency lowering button
5. Indicator
6. Emergency stop
7. Connector for rail roller

Intended use

A ceiling lift provides you with the ability to move a client safely and easily. The lift is used for lifting and moving a client, for example to or from a bed, wheelchair or toilet or from the floor. The apparatus can also be used for walking training, and with the help of some accessories, for balance training and weighing.

HeliQ can be equipped with a Power Traverse that permits motorized horizontal operation and continuous charging direct into the rail.

The ceiling lift has been designed and tested for use indoors and is classified as an IPX4 product.

Important! The product must not be immersed in water.

The climatic conditions should be as follows: ambient temperature from 0 °C to 40 °C, relative ambient humidity from 30% to 80% and ambient pressure from 790 hPa to 1060 hPa.

Product label

The product label contains information about the lift. This label is located on top of the lift.

The product has been tested by an accredited test institute and complies with the necessary requirements in accordance with the Medical Devices Directive for Class 1 products (MDD 93/42/EEC and Regulation 2017/745).

The lift is CE marked and complies with the requirements of EN 10535, EN 60601-1 (Electrical Safety) and EN 60601-1-2 (EMC).

Wireless transmission equipment, mobile telephones, etc. must not be used close to the device as it may affect its function. Particular caution must be observed when using strong sources of interference (such as diathermy, etc.) to ensure that no cables or electrical components are placed on or near the device. If in doubt, ask the person in charge of the equipment or the supplier.

Symbols on the unit and product label

| | |
|--|--|
| | CE marking |
| | Date of manufacture |
| | Manufacturer |
| | Max. load |
| | Instructions: Read the instructions and make sure you understand them before using this product |
| | Article number/reference number "-XX" = country specific code (AU-Australia, GB-Great Britan, EU-Europe, JP-Japan, NA-North America) |
| | Serial number |

| | |
|--|--|
| | Medical Device |
| | Warning |
| | The lift contains lead batteries and electronics and it must be recycled appropriately |
| | Third party certification marking (not applicable for all regions and versions) |
| | Type-B product |
| | Switch the product off/Switch the product on |

Compatibility

Human Care Group slings are designed for use with all Human Care lifts and hanger bars. The products are also compatible with products of other suppliers utilizing the same attachment method for securing slings to a hanger bar and lift. However, due to the wide range of products worldwide, Human Care cannot be responsible for any errors that may result from improper use or application of combinations of slings, hanger bars and lifts of other suppliers. The use of the combined system is at the prescriber's own risk and liability.

To ensure the safety of the client and caregiver, Human Care recommends the following minimum requirements:

1. When using slings with lifts and hanger bars, made by manufacturers other than Human Care, the combined system needs to be designed with the same sling attachment method for the hanger bar of floor lifts and ceiling lifts, meaning loop style slings are combined with a loop style hanger bar, clip style slings are combined with a clip style hanger bar.
2. The involved patient lifts and slings must be CE marked in accordance with Directive 93/42/EEC or Regulation 2017/745 and comply with standard requirements for patient lifts as stated in EN ISO 10535.
3. The caregiver must follow the manufacturer's instruction manuals and recommendations with respect to the use, care, client, caregiver, cleaning, and inspection of the slings and lifts.
4. Different products on the assembled lift system: i.e., lift, hanger bar, sling, scale and other accessories may have different max. loads. It is always the lowest maximum allowable load indicated for the respective products and the system that applies for the entire system. Always check the markings on all individual products in the system.
5. Individual risk assessment, including a practical test of the system and validation with the specific combination used are required to affirm that the size of the sling is correct for the intended use and compatible to the width and design of the hanger bar. The combined system needs to be documented by the competent assessor.

If there are any questions, please contact your local Human Care office or distributor.

Functions

Emergency stop

In an emergency, the emergency stop can be activated by pulling the red belt. To reset the emergency stop, push up the activation plate at the top of the belt.

Emergency lowering

In the event of an emergency, the emergency lowering button on the lift or hand control can be pressed to quickly lower the lift. An audio signal indicates that the lift is lowering. Always ensure that emergency lowering can be performed as safely as possible for the client. The

emergency lowering button only works when the lift is switched on and emergency stop has been deactivated. When the emergency lowering is active, the end position switches are inactive. If the button is pushed until the lift is in its lowest position, the lift will start to go up again. Always stop at the lowest position.

The emergency lowering button must only be used in emergency situations.

Safety catch

The lift is equipped with a safety catch. This prevents involuntary falls. If the lift belt is pulled out too quickly, the safety catch will engage, thus preventing the fall.

Mechanical emergency lowering

During power loss and the hoist is loaded, user in the hoist, the mechanical emergency

lowering is activated and the hoist is lowering in a safe and reliable manner.

Overheating

The lift has overheating protection which stops the motor if it gets too hot. Overheating can occur if the motor is overloaded or operated continuously for quite a long time. The lift's work cycle (work:pause) is 15:85. Maximum continuous operation is two minutes. If this is exceeded, there is a risk of the motor overheating.

Belt tension guard


The lift is equipped with a belt tension guard that prevents unintentional discharging of the lift belt. It must be loaded in order to be able to pull the belt out. Pull gently on the belt while at the same time pressing the Down button on the hand control or control panel. The weight of the bar is sufficient for the belt to be fed out.

Suspending the lift from the rail

HeliQ is mounted directly in the rail roller, which is located in the rail. Do not use any distance strap or similar.

Remove the pins (A) holding the two bolts (B) to the rail roller. Remove the bolts and raise the lift to the rail roller. Make sure that the holes in the connector on the lift and the holes in the rail roller are lined up, and replace the bolts in the holes. Secure the bolts with the two pins. The lift can be mounted in either direction. Switch on at the mains and the lift is ready for use.

Below image: Correct placement of the locking pins.

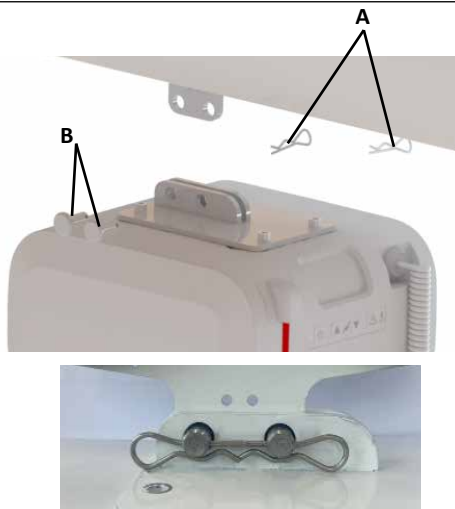Hand control

The hand control is equipped with an Up button (2) and a Down button (3) for raising and lowering the lift. The hand control is also equipped with an OFF/ON button (1) for activation/deactivation of the lift and an emergency lowering button (4) for quickly lowering the lift.


Push the emergency lowering button together with the down button at least for 3 seconds simultaneously to activate the emergency lowering. Older lifts might not have the dual control and the delay.

Lifts equipped with a Power Traverse have a hand control with seven buttons. In addition to the buttons described above, there are also buttons for moving the lift forwards (7) and backwards (6), as well as an alarm button (5).

When the lift is loaded with less than approx. 40 kg it can be run at two speeds. If the Up or Down button is pressed for longer than 3 seconds, the lift will run in the faster mode. If

the lift is loaded with more than approx. 40 kg it will only run in the slower mode.

On the short side of the hand control there is an input for charging the lift and an input (USB-C) for communicating with the lift.

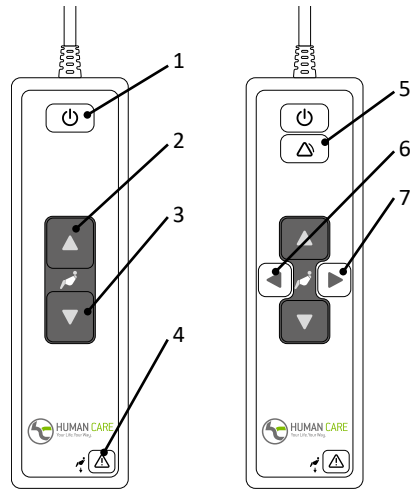HeliQ


HeliQ Power traverse

Secondary control panel

The lift is equipped with a secondary control panel located on the lift itself. This control panel can be used as an alternative to the hand control. The secondary control panel is equipped with an Up button (2) and a Down button (3) for raising and lowering the lift. The control panel is also equipped with an OFF/ON button (1) for activation/deactivation of the lift and an emergency lowering button (4) for quickly lowering the lift.

Press emergency lowering button and down button at least 3 seconds simultaneously to activate the emergency lowering. Older lifters may not have double key command and delay.

When the lift is loaded with less than approx. 40 kg it can be run at two speeds. If the Up or Down button is pressed for longer than 3 seconds, the lift will run in the faster mode. If the lift is loaded with more than approx. 40 kg it will only run in the slower mode. The force required to operate the push buttons is 4.4N.

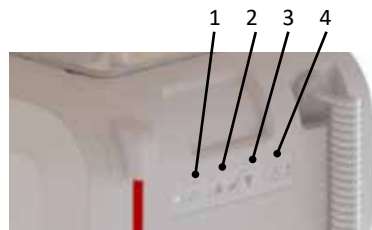Indicator

The indicator light has the following indications, which are only accessible when the lift is switched on.

| | |
|-----------------|---|
| Flashing green | The lift is on and ready for use. Battery level is Full. |
| Flashing yellow | The lift is on and ready for use. Battery level is Normal. It might be a good idea to charge up the lift. |
| Flashing red | The lift is on and ready for use. Battery level is Low. The lift must be charged up immediately. |
| Constant green | The lift is in operation or loaded. |
| Constant red | The lift safety catch has engaged. (See troubleshooting table). |

| | |
|-----------------------|---|
| Constant yellow | The lift is charging. |
| Flashing green/yellow | The lift's service interval has passed. Service and annual inspection must be performed immediately. |
| Flashing red/yellow | The lift's service interval has passed. Service and annual inspection must be performed immediately. Battery level is Low. The lift must be charged up immediately. |

Charger

The indicator on the lift flashes red when the batteries need recharging. Plug the charger cable into the hand control or put the hand control into the charger cradle and charging will start. When the lift is charging the indicator glows yellow.

The lift can only be used with a charger recommended by Human Care, article number 50880.

- Charge the lift regularly, preferably every night
- The charger must be disconnected when you use the lift
- The unit charges regardless of the position of the main switch
- If the emergency stop is activated the lift cannot be charged. Check that the emergency stop is not activated when charging.
- The charger must only be connected to an earthed socket.

If rail charging is available via the rail, the lift will charge continuously when the lift is not in use and is within the charging zone. When rail charging is available, the lift does not need to be charged via the hand control. Get in touch with the system owner or Human Care if you are uncertain as to what type of charging to use.

Important! When the batteries reach a critically low battery level, it will no longer be possible to perform more lifts (up or down). However, there is still enough battery power for emergency lowering. See separate manual for the charger, 99599.

Press the Up button on the hand control or control panel on the lift to lift the client. Press the Down button on the hand control or control panel to lower the client. Pressing the button activates a soft start, and releasing it causes a soft stop.

Detaching the lift from the rail system

Remove the pins holding the two bolts to the rail roller. Hold the lift firmly and remove both bolts. The lift is now disconnected from the rail roller.

Place the lift on a flat surface. Press the Up button on the lift to move the belt back. Make sure

the belt is kept taut and gets rolled up in the lift straight, without folding.

Important! When the lift is taken down, the belt must be kept taut while it is fed into the lift.

Using the lift

Human Care offers several different types of bar in a variety of widths. The bar chosen should be appropriate for the client's needs and the type of sling used. Place the client in the prescribed

sling according to the sling instructions and hang the sling in the hooks on the bar. Check that the sling's lifting straps are properly anchored in all hooks.

Data communication

HeliQ is able to communicate with a standard PC (Win 7/Win 8). This communication makes it possible to obtain statistics and other information on the lift.

To communicate with the ceiling lift you need the Human Care Data and Service Application, software that can be downloaded from the Human Care website, www.humancaregroup.com.

Connect the lift to your PC with a standard USB-C cable. Start the Human Care Data and Service Application and follow the instructions. The lift must be switched on for communication to be possible. The lift is deactivated while communication takes place.

The following information can be read from the lift.

- Article number, serial number and production date
- Max. user weight
- Software version
- Service dates
- Periodic usage statistics
- Error log

Connected computer equipment must fulfil the requirements of IEC 60950-1 or IEC 60601-1

The latest software for the lift should always be used. You can find the latest version at our website in the product section.

Hanger Bars

Connect the hanger bar to the hook on the lift belt by threading the hook diagonally from beneath (see fig).

- 72754 Hanger bar HeliQ 2P 35 cm
- 72755 Hanger bar HeliQ 2P 45 cm
- 72757 Hanger bar HeliQ 2P 60 cm
- 92938 Hanger bar HeliQ 4P

Lifting height

The lifting height corresponds to the difference between the highest and lowest position the lift can reach. This corresponds to the length of the lift belt. The length of the lift belt is 2120mm (83,5")/1720 mm (66.1").

If you want to extend the lift belt there are distance straps in various lengths available as extras. The lifting height of the lift is not affected by the length of the distance strap. The total height of the system consists of the rail in which the lift is mounted, the rail roller in the rail, a distance strap if relevant and the lift.

Accessories

If the maximum weight for any of the accessories, e.g. slings, lift bars and any other accessories is NOT the same as for the lift, then the lowest max. weight given for any one product is the one that shall apply. Always check the label on the lift and lifting accessories or contact Human Care if you have any questions or if anything is unclear.

Slings

Human Care's range of slings is wide enough to meet the user's various requirements and needs. Slings are tested in accordance with Human Care's test log, which can be downloaded from www.humancaregroup.com. A test must be performed for each client individually and in accordance with their needs. Consult the user manual for the relevant sling to make sure the product is used correctly.

Human Care recommends using only Human Care slings with Human Care lifts. Visit www.humancaregroup.com for further information about which models and sizes are offered.

Motorized operation

HeliQ offers the ability to connect to motorized operation. Motorized operation allows the lift to be moved along the rail with minimal effort from the caregiver. Motorized operation needs to be installed by a technician with Human Care qualifications.

The power traverse is equipped with overheating protection which stops the motor if it overheats. Overheating can occur if the motor is running continuously during a longer period. Continuous operation is maximum 5 minutes. More than this will lead to overheating.

Troubleshooting

| Description | Indicator | Action |
|---|---|--|
| Nothing happens when you press the Up button or the Down button on the hand control or the control panel | Not lit | <ol style="list-style-type: none"> 1) Check that the mains switch is on. 2) The batteries are fully discharged. Charge the lift immediately 3) The emergency stop is pushed in. Deactivate by turning the red button clockwise. |
| Nothing happens when you press the Up button or the Down button on the hand control | Flashing green Flashing yellow Flashing red | Check the hand control. If the hand control is defective it needs to be replaced. Contact a service technician. |
| Nothing happens when you press the Up button or the Down button on the hand control or the control panel. | Constant red | <ol style="list-style-type: none"> 1) Lift motor is overheated. Wait a while and try again. 2) Lift is overloaded. Lower using the emergency lowering button |
| Nothing happens when you press the Up button or the Down button on the hand control or the control panel. | Constant yellow | Charging in progress. Interrupt charging. |
| Nothing happens when you press the Up button on the hand control or the control panel. | Flashing green Flashing yellow Flashing red | The lift is in its maximum upper position. Press the Down button on the hand control or the control panel to lower the lift. |
| Nothing happens when you press the Down button on the hand control or the control panel. | Flashing green Flashing yellow Flashing red | <ol style="list-style-type: none"> 1) The lift is in its maximum lower position. Press the Up button on the hand control or the control panel to raise the lift. 2) Lift is unloaded. Pull gently on the hook on the lift belt or hang up a lift bar and press the Down button on the hand control or control panel to lower it. |
| The indicator switches from flashing green or yellow to flashing red. | Flashing red | The batteries need charging but you can still lift the client a few times. The lift should be set to charge as soon as possible. |

Care instructions for the lift and accessories

Care instructions for the lift

- The lift should be recharged every night
- Wipe the ceiling lift with a slightly dampened cloth with a water solution using a mild detergent.
- Never immerse the lift in water, or use any other fluids to clean the lift.
- NEVER use strong detergents on the ceiling lift.
- If the lift needs to be disinfected, Du-Pont™ Virkon® can be used
- The lift belt must be washed and disinfected as instructed above
- Activate the emergency stop if the ceiling lift is not to be used for a few days, otherwise the batteries will run down
- The lift must always be stored or moved lying horizontally

Annual inspection

The lift must be inspected at least once a year by someone with Human Care qualifications and in accordance with Human Care's directive. Pay particular attention to parts that show wear. The lift's service indicator will flash yellow when a service needs to be performed.

Daily inspection

- Always check the sling before use. There must be no loose seams or other signs of wear to lift belts, belts and fabric. If there are, replace the sling
- Check the ceiling lift for visible signs of damage
- Check that the lift is correctly mounted in the distance strap or rail roller

Monthly inspection

Check the lift belts for damage or wear and tear. Feed out enough of the lift belts to allow their

entire length to be checked. Contact a retailer if the belt is damaged.

Service

Installation, service and maintenance must only be performed by personnel with Human Care qualifications and in accordance with Human Care's directive. Only original spare parts may be used.

Service agreements

Human Care invites you to sign a service agreement for annual inspection and test loading.

Transport and storage

During transport, or when the lift is not going to be used for some time, the mains switch must be off or the emergency stop button pushed in, otherwise the batteries will self-discharge. The lift must be transported and stored lying horizontally. The climatic conditions should be as follows: ambient temperature from 0 °C to 40 °C, relative ambient humidity from 30% to 80% and ambient pressure from 790 hPa to 1060 hPa.

Product lifetime

The product's expected lifetime is ten (10) years or 10.000 lifts, if the product is used as intended and maintained according to the manufacturer's instructions, depending on the intensity of use and maximum load applied during use. If the product label is no longer legible, the product should be discarded.

Warranty and support

If you need information or support, please contact your Human Care office, www.humancare-group.com, or distributor.

Guidance and manufacturer's declaration

| Guidance and manufacturer's declaration – electromagnetic emissions | | | |
|---|---|---|---|
| The Human Care lift is intended for use in the electromagnetic environment specified below. The customer or the user of the Human Care lift should assure that it is used in such an environment. | | | |
| Emission test | Compliance | Electromagnetic environment - guidance | |
| RF emissions CISPR 11 | Group 1 | The Human Care lift uses RF energy only for its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment. | |
| RF emissions CISPR 11 | Class B | The Human Care lift is suitable for use in all establishments including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes. | |
| Harmonic emissions IEC 61000-3-2 | Class A | | |
| Voltage fluctuations/flicker emissions IEC 61000-3-3 | Complies | | |
| Guidance and manufacturer's declaration – electromagnetic immunity | | | |
| The Human Care lift is intended for use in the electromagnetic environment specified below. The customer or the user of the Human Care lift should assure that it is used in such an environment. | | | |
| Immunity test | IEC 60601 test level | Compliance level | Electromagnetic environment - guidance |
| Electrostatic discharge (ESD) IEC 61000-4-2 | +/- 6 kV contact +/- 8 kV air | +/- 6 kV contact +/- 8 kV air | Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %. |
| Electrical fast transient / Burst IEC 61000-4-4 | +/- 2 kV for power supply lines +/- 1 kV for input/output lines | +/- 2 kV for power supply lines n/a. for input/output lines | Mains power quality should be that of a typical commercial or hospital environment. |
| Surge IEC 61000-4-5 | +/- 1 kV differential mode +/- 2 kV common mode | +/- 1 kV differential mode n/a. for common mode | Mains power quality should be that of a typical commercial or hospital environment. |
| Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11 | <5 % UT (>95 % dip in UT) for 0,5 cycle 40 % UT (60 % dip in UT) for 5 cycles 70 % UT (30 % dip in UT) for 25 cycles <5 % UT (>95 % dip in UT) for 5 sec | <5 % UT (>95 % dip in UT) for 0,5 cycle 40 % UT (60 % dip in UT) for 5 cycles 70 % UT (30 % dip in UT) for 25 cycles <5 % UT (>95 % dip in UT) for 5 sec | Mains power quality should be that of a typical commercial or hospital environment. If the user of the Human Care lift requires continued operation during power mains interruptions, it is recommended that the lift be powered from an uninterruptible power supply or battery. |

| | | | |
|--|-------|-------|--|
| Power frequency (50/60 Hz) magnetic field IEC 61000-4-8 | 3 A/m | 3 A/m | Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment |
|--|-------|-------|--|

NOTE UT is the a.c. mains voltage prior to application of the test level.

Guidance and manufacturer's declaration – electromagnetic immunity

The Human Care lift is intended for use in the electromagnetic environment specified below. The customer or the user of the Human Care lift should assure that it is used in such an environment.

| Immunity test | IEC 60601 test level | Compliance level | Electromagnetic environment- guidance |
|-------------------------------|-----------------------------|------------------|--|
| Conducted RF IEC 61000-4-6 | 3 Vrms 150 kHz to 80 MHz | 3 Vrms | <p>Portable and mobile RF communications equipment should be used no closer to any part of the Human Care lifts, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.</p> <p>Recommended separation distance $d = [3,5/3] \sqrt{P}$</p> <p>$d = [3,5/3] \sqrt{P}$ 80 MHz to 800 MHz</p> <p>$d = [7/3] \sqrt{P}$ 800 MHz to 2,5 GHz</p> <p>where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m).</p> <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, should be less than the compliance level in each frequency range. b</p> <p>Interference may occur in the vicinity of equipment marked with the following symbol.</p>  |
| Radiated RF IEC 61000-4-3 | 3 V/m 80MHz to 2,5GHz | 3 V/m | |
| | 10 V/m 800MHz to 2,5GHz | 10V/m | |

NOTE 1 At 80MHz and 800MHz, the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflected from structures, objects and people.

a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the Human Care lifts is are used exceeds the applicable RF compliance level above, the Human Care lifts should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the Human Care lifts.

b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 10 V/m.

Recommended separation distances between portable and mobile RF communications equipment and the Human Care lifts.

The Human Care lifts are intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the Human Care lifts can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the Human Care lifts as recommended below, according to the maximum output power of the communications equipment.

| Rated maximum output power of transmitter W | Separation distance according to frequency of transmitter m | | |
|--|--|--------------------------------|---------------------------------|
| | 150 kHz to 80 MHz d = 1,2√P | 80 MHz to 800 MHz d = 1,2√P | 800 MHz to 2.5 GHz d = 0,7√P |
| 0.01 | 0.12 | 0.12 | 0.07 |
| 0.1 | 0,37 | 0,37 | 0,22 |
| 1 | 1,16 | 1,16 | 0,7 |
| 10 | 3,67 | 3,67 | 2,21 |
| 100 | 11.6 | 11.6 | 7 |

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

Note 1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

Recycling

Once the product has reached the end of its life, it should be disposed according to your country's rules, especially for the separate collection of electrical and electronic products.

Correct disposal helps prevent negative consequences for the environment and human health.

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