



Now eligible for ADP funding!

AA0001810



Smart and slim design

neXus oNe, is one of our rollators with the slimmest design, making it perfect in narrow spaces.



Easy to fold

The best selling neXus series of rollators have innovated the way we think about walking aids. With the first-ever patented cross-folding frame, it set the standard for future walking aids design. Like the neXus 3 rollator, the neXus oNe rollator features a cross-folding frame, making it exceptionally easy to fold. The soft zippered bag folds together while remaining on the rollator. It can stand independently in a folded position with the bag on the rollator. Although the neXus oNe rollator is primarily for indoor use, its 6" wheels makes it a great rollator to navigate in smaller spaces such as at home, in the aisles of a store or an airport. The cross-folding design also enables the rollator to be partially folded to pass through narrow spaces.

Cable-free and anatomic brake


The cable-free and anatomic brake handles increases safety and comfortability. No cables catch on objects around the home. The cable-free braking system also ensures years of trouble-free use. We have had anatomy and strong visual communication in mind when designing the brake handles. They are comfortable, and the red ball makes the brake handles visible. The cross-folding structure enables superior walking ergonomics with ample walking space. The padded, contoured, and cushioned seat offers a comfortable resting spot.

Available in three sizes

The neXus oNe is available in three different seat heights; neXus oNe 17, neXus oNe 20, and neXus oNe 22.



Features

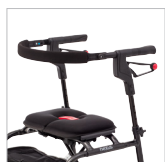


Padded Backrest

Extra back comfort while seated. Can be flipped both ways and is included with the rollator. Part no. P48010.



Foam Cushioned Seat

Offers a comfortable resting spot while walking.

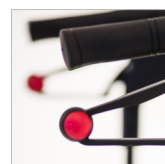

Anatomic Brake Handles

The anatomic brake handles are easy and comfortable to use.

Accessories




Tray

Practical for resting lightweight items. Maximum load 2 kg/4.4lb. Part no. 93316.



Hard Basket

Basket ideal for shopping. Maximum load 4,5kg/10 lb. Part no. 4510.

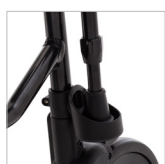

Cane Holder

Simplifies transport of a cane or umbrella. Part no. 93317.



Cup Holder

Attaches in seconds without tools. Fits standard size bottle. Part no. 90394.



Curb Climper

Allows the rollator to be slightly lifted onto a curb. Part no. 93315.



Identification Tag Green

Attaches to the rollator in seconds. Part no. 90393.



Identification Tag Canada

Attaches to the rollator in seconds. Part no. 90392.



Soft zippered bag

Soft bag that folds with the mobility aid. Maximum capacity 5kg/11lb. Part no. 93280.

Technical specifications

Part Name	Part no.	Handle Height (Min - Max)	Seat Height	Overall Width	Width between Handles	Distance between Rear Wheels	Wheel Size	Product Weight
neXus oNe 17	490017	72 cm - 80 cm 28.4" - 31.5"	47 cm 17"	58 cm 23"	44 cm 17.3"	44 cm 17.3"	15 cm x 3 cm 6" x 1.2"	7.7 kg 17lb
neXus oNe 20	490020	74 cm - 89 cm 29" - 35"	51 cm 20"	58 cm 23"	44 cm 17.3"	44 cm 17.3"	15 cm x 3 cm 6" x 1.2"	7.9 kg 17.4lb
neXus oNe 22	490022	79 cm - 97 cm 31" - 38"	56 cm 22"	58 cm 23"	44 cm 17.3"	44 cm 17.3"	15 cm x 3 cm 6" x 1.2"	8.1 kg 17.9lb

Maximum User Weight: 160kg / 353lb

” We improve freedom of movement and quality of care for people.

