



Cableless

Comfortable & Stylish

High Weight Capacity!



Our bestselling mobility aid

neXus 3, is our bestselling rollator with a wide backstrap and robust wheels, making it a perfect everyday-life-companion.

Easy to fold

With the first-ever patented cross-folding frame, it set the standard for future walking aids design. The neXus 3 rollator features a cross-folding frame, making it exceptionally easy to fold. The soft zippered bag folds together while remaining on the rollator. It can stand independently in a folded position with the bag on the rollator. The cross-folding design also enables the rollator to be partially folded to pass through narrow spaces.

Cableless and wide backstrap

The cable-free and anatomic brake handles increases safety and comfort. No cables catch on objects around the home. The cable-free braking system also ensures years of trouble-free use. We have had anatomy and strong

visual communication in mind when designing the brake handles. They are comfortable, and the red ball is designed for those with limited hand dexterity and allows the user to apply/disengage the brakes, using the palm of their hands. The cross-folding structure enables superior walking ergonomics with ample walking space.

The padded, contoured, and cushioned seat offers a comfortable resting spot. The backstrap can be flipped to allow for two-position seating as well as a visual aid for achieving optimal gait.

Available in black, blue and red

The neXus 3 is available in three different sizes; neXus 3 Super Low (19"), neXus 3 Low (22"), and neXus 3 Standard (24")*. All sizes are available in Black, Blue and Red (Burgundy) colour. The neXus 3 Wide are only available in Black.

* Measurement for sizes is floor to bottom of seat pan.



Mobility Aids



Healthcare Beds



Lifting Solutions




Hygiene



Convertible Chairs

Features

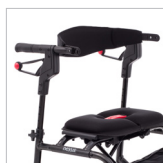


Padded Backrest

Back comfort while seated. Can be flipped both ways and is included with the rollator.



Foam Cushioned Seat

Offers a comfortable resting spot while walking.

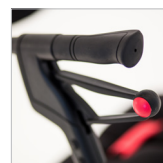

Anatomic Brake Handles

The anatomic brake handles are easy and comfortable to use.

Accessories



Hard Basket

Basket ideal for shopping. Maximum load 4.5kg/10 lb. Part no. 4510.



Tray

Practical for resting lightweight items. Maximum load 2 kg/4.4lb. Part no. 93316.




Curb Climber

Allows the rollator to be slightly lifted onto a curb. Part no. 93315.




Cane or Umbrella Holder

Simplifies transport of a cane or umbrella. Part no. 93317.

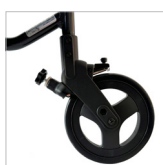

Slow Down Brake

Enables controlled walking speed. Delivered in pairs. Part no. 93122.

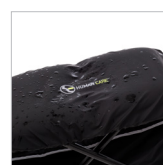

Seat Cover

Waterproof seat cover that is easy to attach to the seat. Part no. 93313.



Cup Holder

Attaches in seconds without tools. Fits standard size bottle and cups. Part no. 90394.



Identification Tag Canada

Attaches to the rollator in seconds. Part no. 90392.



Soft Zippered Bag

Soft bag that folds with the mobility aid. Maximum capacity 5kg /11lb. Part no. 93312.

Technical specifications

Part Name	Part no.	Handle Height (Min - Max)	Seat Height	Overall Width	Width between Handles	Distance between Rear Wheels	Wheel Size	Product Weight	Maximum user Weight
neXus 3 Super Low Black Blue Red (Burgundy)	470018* 470018B 470018R	78cm - 86cm / 30.5" - 34.1"	48cm / 19"	58cm / 23"	45cm / 17.5"	47cm / 18.5"	20cm x 3cm / 8" x 1.2"	8.39 kg / 18.5lb	200 kg / 441 lb
neXus 3 Low Black Blue Red (Burgundy)	470022 470022B 470022R	78cm - 93cm / 30.5" - 36.5"	56cm / 22"	58cm / 23"	45cm / 17.5"	47cm / 18.5"	20cm x 3cm / 8" x 1.2"	8.48 kg / 18.7lb	200 kg / 441 lb
neXus 3 Standard Black Blue Red (Burgundy)	470024 470024B 470024R	83cm - 101cm / 32.5" - 39.8"	61cm / 24"	58cm / 23"	45cm / 17.5"	47cm / 18.5"	20cm x 3cm / 8" x 1.2"	8.62 kg / 19.0lb	200 kg / 441 lb
neXus 3 Low Wide Black	479022	78cm - 93cm / 30.5 - 36.5"	56cm / 22"	67cm / 26.4"	53cm / 20.9"	56.5cm / 22.2"	20cm x 3cm / 8 x 1.2"	8.4 kg / 18.5lb	250kg / 551lbs
neXus 3 Standard Wide Black	479024	83cm - 101cm / 32.5 - 39.8"	61cm / 24"	67cm / 26.4"	53cm / 20.9"	56.5cm / 22.2"	20cm x 3cm / 8 x 1.2"	9 kg / 19.8lb	250kg / 551lbs

* Article number 470018S (neXus 3 Super Low Black) can be equipped with a different handle kit to lower the handle height to 27,5" - 30,5".



*We improve freedom of movement
and quality of care for people.*

Human Care Canada
10-155 Colonnade Road
Ottawa, ON K2E 7K1, CANADA
+1 613 723 6734
info.ca@humancaregroup.com

www.humancaregroup.com