



Comfortable

Narrow Design

Head Support



Silhouette Sling Soft
Art. Nr: 25019

The Silhouette Sling is a comfortable sling without reinforcements, handles or padding. It adapts to the patient's body to offer good support for the entire body. The sling can be used in most lifting situations and is especially suitable for transferring to and from wheelchairs, molded seats or seating shells.



Silhouette Sling Polyester Net
Art. Nr: 25015

Two Fabric Choices

The Silhouette Sling is available in either Polyester Net or Soft Air Mesh, both of which allow moisture and heat to transfer efficiently from the patient. These slings can be left behind the patient for a short period of time after the transfer is complete. The Polyester Net fabric is flexible, airy and dries quickly which makes it suitable for use in bathing environments.

s available in sizes Junior Small to Extra large, accommodating a wide assortment of users.

Easy Care and Easy Fit

Our slings are manufactured to the highest quality standards. They feature reinforced straps and high strength nylon webbing around the sling to make them strong enough to resist even the most active use. Each sling label is colour coded to ensure quick and easy selection of the correct sling and the straps are colour coded for easy determination of required positioning.

Narrow Overall Width

The Silhouette Sling has a narrow overall width which makes it an ideal option for users who are smaller in stature. The Silhouette Sling can be used for most lifting situations with a 2-point, 3-point or 4-point suspension and is suitable for patients with poor trunk stability. The Silhouette Sling



Mobility Aids



Healthcare Beds



Lifting solutions



Hygiene



Convertible Chairs

Accessories



Extension Loops
Used to extend lift straps in order to alter lifting position or to hold leg parts together. Assorted lengths. Part no.:21180-*



Leg Padding and Supports
Ideal for users with sensitive skin or pain. Part no.: 21146-S, M, L Reinforced supports. One size. Part no. 25102



Upholstered Cover
Protects user from hanger bar impact. Not for use in a wet environment. Part no. 50866



Hanger Bars for 2-, 3- or 4-point Suspension
Part no. 57551, 55764 (2- or3-Point - Quick Release), 55619 (4-Point - Attach to leg hook)


Lifting Situation



Use to/from a bed
Suitable for transfers and repositioning to and from a bed. Can be left behind patient for short period of time.




Use to/from a chair
Suitable for transfers and repositioning to and from a chair. Can be left behind patient for short period of time.

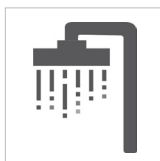

Use for showering
Only Polyester Net fabric can be used for this application.

Care and Features

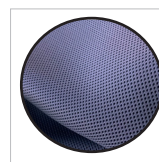

Soft Air Mesh
Lightweight, cool and breathable material. Suitable for sensitive skin.



Polyest Net
Soft, thin material that breathes and dries quickly. Ideal for bathing or showering.




Care
Care instructions located on sling label. Use a hospital grade laundry bag while laundering.

Technical specifications

| Article number | | | Size | Colour | A | B | C |
|-------------------------|--------------------------|---------------|---------------|--------|---------------|---------------|----------------|
| Polyester Net | Polyester Net padded leg | Soft Air Mesh | | | | | |
| 25015-JS | 25018-JS | 25019-JS | Junior Small | JS | 18 cm / 7,1" | 62 cm / 24,4" | 64 cm / 25,2" |
| 25015-JM | 25018-JM | 25019-JM | Junior Medium | JM | 22 cm / 8,7" | 67 cm / 26,4" | 71 cm / 28" |
| 25015-S | 25018-S | 25019-S | Small | S | 26 cm / 10,2" | 77 cm / 30,3" | 79 cm / 31,1" |
| 25015-M | 25018-M | 25019-M | Medium | M | 29 cm / 11,4" | 87 cm / 34,2" | 90 cm / 35,4" |
| 25015-L | 25018-L | 25019-L | Large | L | 35 cm / 13,8" | 93 cm / 36,6" | 96 cm / 37,8" |
| 25015-XL | 25018-XL | 25019-XL | Extra Large | XL | 40 cm / 15,8" | 97 cm / 38,2" | 106 cm / 41,8" |
| maximales Nutzergewicht | | | 300kg / 660lb | | | | |

* Available in assorted lengths: 10cm/3,9" (21180-10), 15cm/5,9" (21180-15), 20cm/7,8" (21180-20), 25cm/9,8" (21180-25), 35cm/13,7" (21180-35)

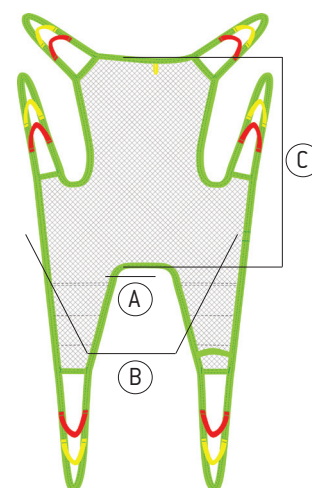

99020 Rev.3, EN, Product Sheet, Silhouette Sling, February 2022



We improve freedom of movement and quality of care for people.

www.humancaregroup.com



HUMAN CARE SWEDEN (HQ)
Årstaängsvägen 21B
117 43 Stockholm
+46 8 665 35 00
info@humancaregroup.com

HUMAN CARE GERMANY
Industriestraße 13, Gebäude 1/9
63755 Alzenau
+49 (0) 6023 994 936 8
infode@humancaregroup.com

HUMAN CARE CANADA
10-155 Colonnade Road
ON K2E 7K1 Ottawa
+1 613 723 6734
info.ca@humancaregroup.com

HUMAN CARE THE NETHERLANDS
Elspeterweg 124
8076 PA Vierhouten
+31 577 41 21 71
info.nl@humancaregroup.com

HUMAN CARE UNITED STATES
8006 Cameron rd, Suite K
TX 78754 Austin
+1 5 124 76 7199
info.us@humancaregroup.com