

# Silhouette Sling



Comfortable

Narrow Design

Head Support



The Silhouette Sling is a comfortable sling without reinforcements, handles or padding. It adapts to the patient's body to offer good support for the entire body. The sling can be used in most lifting situations and is especially suitable for transferring to and from wheelchairs, molded seats or seating shells.



## Two Fabric Choices

The Silhouette Sling is available in either Polyester Net or Soft Air Mesh, both of which allow moisture and heat to transfer efficiently from the patient. These slings can be left behind the patient for a short period of time after the transfer is complete. The Polyester Net fabric is flexible, airy and dries quickly which makes it suitable for use in bathing environments.

## Easy Care and Easy Fit

Our slings are manufactured to the highest quality standards. They feature reinforced straps and high strength nylon webbing around the sling to make them strong enough to resist even the most active use. Each sling label is colour coded to ensure quick and easy selection of the correct sling and the straps are colour coded for easy determination of required positioning.



## Narrow Overall Width

The Silhouette Sling has a narrow overall width which makes it an ideal option for users who are smaller in stature. The Silhouette Sling can be used for most lifting situations with a 2-point, 3-point or 4-point suspension and is suitable for patients with poor trunk stability. The Silhouette Slings are available in sizes Junior Small to Extra large - accommodating a wide assortment of users with different needs.



Mobility Aids



Healthcare Beds



Lifting solutions



Hygiene



Convertible Chairs

## Accessories



### Extension Loops

Used to extend lift straps in order to alter lifting position or to hold leg parts together. Assorted lengths. Part no.:21180-\*



### Leg Padding and Supports

Ideal for users with sensitive skin or pain. Part no.: 21146-S, M, L Reinforced supports. One size. Part no. 25102



### Upholstered Cover

Protects user from hanger bar impact. Not for use in a wet environment. Part no. 50866



### Hanger Bars for 2-, 3- or 4-point Suspension

Part no. 57551, 55764 (2- or 3-Point - Quick Release), 55619 (4-Point - Attach to leg hook)

## Lifting Situation




### Use to/from a bed

Suitable for transfers and repositioning to and from a bed. Can be left behind patient for short period of time.



### Use to/from a chair


Suitable for transfers and repositioning to and from a chair. Can be left behind patient for short period of time.

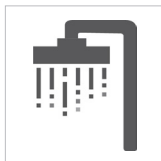

### Use for showering

Only Polyester Net fabric can be used for this application.

## Care and Features

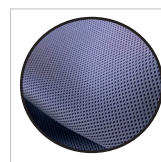

### Soft Air Mesh

Lightweight, cool and breathable material. Suitable for sensitive skin.



### Polyest Net

Soft, thin material that breathes and dries quickly. Ideal for bathing or showering.




### Care

Care instructions located on sling label. Use a hospital grade laundry bag while laundering.

## Technical specifications

Article number			Size	Colour	A	B	C
Polyester Net	Polyester Net padded leg	Soft Air Mesh					
25015-JS	25018-JS	25019-JS	Junior Small	JS	18 cm / 7,1"	62 cm / 24,4"	64 cm / 25,2"
25015-JM	25018-JM	25019-JM	Junior Medium	JM	22 cm / 8,7"	67 cm / 26,4"	71 cm / 28"
25015-S	25018-S	25019-S	Small	S	26 cm / 10,2"	77 cm / 30,3"	79 cm / 31,1"
25015-M	25018-M	25019-M	Medium	M	29 cm / 11,4"	87 cm / 34,2"	90 cm / 35,4"
25015-L	25018-L	25019-L	Large	L	35 cm / 13,8"	93 cm / 36,6"	96 cm / 37,8"
25015-XL	25018-XL	25019-XL	Extra Large	XL	40 cm / 15,8"	97 cm / 38,2"	106 cm / 41,8"
Maximal user weight			300kg / 660lb				

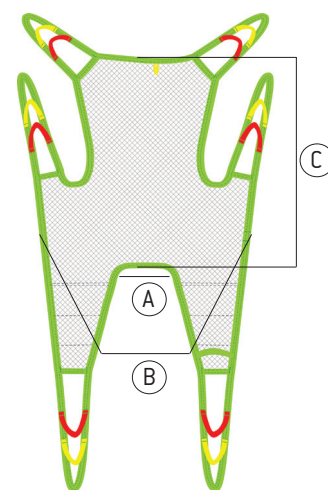

\* Available in assorted lengths: 10cm/3,9" (21180-10), 15cm/5,9" (21180-15), 20cm/7,8" (21180-20), 25cm/9,8" (21180-25), 35cm/13,7" (21180-35)



*We improve freedom of movement  
and quality of care for people.*

[www.humancaregroup.com](http://www.humancaregroup.com)



### HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B  
117 43 Stockholm  
+46 8 665 35 00  
info@humancaregroup.com

### HUMAN CARE CANADA

10-155 Colonnade Road  
ON K2E 7K1 Ottawa  
+1 613 723 6734  
info.ca@humancaregroup.com

### HUMAN CARE THE NETHERLANDS

Elsperweg 124  
8076 PA Vierhouten  
+31 577 41 21 71  
info.nl@humancaregroup.com

### HUMAN CARE UNITED STATES

8006 Cameron rd, Suite K  
TX 78754 Austin  
+1 512 76 7199  
info.us@humancaregroup.com