

highback plusline sling



The Highback Plusline Sling is a high back model with divided leg parts. It has a generous overall width to offer a stretched seating position that is well-suited for heavier patients. The reclined seating offers a balanced weight distribution and the combined Polyester/Polyester Net fabric makes it easy to use and apply.

Comfortable Seating

The Highback Plusline features reinforced and divided pressure relieving leg parts as well as a reinforced head support. Combined with the generous overall width of this sling, it makes the Highback Plusline an extremely comfortable sling to use. The reclined position evenly distributes weight throughout the sling ensuring that the user does not experience any discomfort during the lifting process.

Bariatric Capacity

The Highback Plusline is designed to lift patients of up to 318kg/701lb. Although this sling can be used with both 2-Point and 3-Point suspension, it is recommended that

4-Point suspension is used as it provides a stretched seating position that is well-suited for heavier patients.

Easy Care and Easy Fit

Our slings are manufactured to the highest quality standards. They feature reinforced straps and high strength nylon webbing around the sling to make them strong enough to resist even the most active use. Each sling label is colour coded to ensure quick and easy selection of the correct sling and the straps are colour coded for easy determination of required positioning.

HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B
117 43 Stockholm
Phone: +46 8 665 35 00
Fax: +46 8 665 35 10
info.se@humancaregroup.com

HUMAN CARE AUSTRALIA

Unit 4, 49 Bertie Street
Port Melbourne, VIC 3207
Phone : +61 1300 132 736
Fax : + 61 (03) 9645 4954
info.au@humancaregroup.com

HUMAN CARE CANADA

10-155 Colonnade Road
Ottawa, ON K2E 7K1
Phone: +613.723.6734
Fax: +613.723.1058
info.ca@humancaregroup.com

HUMAN CARE NEDERLAND

Elspeterweg 124
8076 PA Vierhouten
Phone: +31 577 412 171
Fax: +31 577 412 170
info.nl@humancaregroup.com

HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K
Austin, TX 78754
Phone: +512.476.7199
Fax: +512.476.7190
info.us@humancaregroup.com

Accessories



Extension Loops

Used to extend lift straps in order to alter lifting position or to hold leg parts together. Assorted lengths. Part no. 21180.*



Upholstered Cover

Protects user from hanger bar impact. Not for use in a wet environment. Part no. 50866



4-Point Suspension

Enables 4-point lifting. Part no. 70800



Hanger Bars for 2-, 3- or 4-point Suspension

Part no. 57551, 55764 [2- or 3-Point - Quick Release], 55619 [4-Point - Attach to leg hook]

Lifting Situations



Use to/from a bed

Suitable for transfers and repositioning to and from a bed. Can be left behind patient for short period of time.



Use to/from a chair


Suitable for transfers and repositioning to and from a chair. Can be left behind patient for short period of time.



Use to/from the floor


Suitable for transfers from the floor.

Care and Features

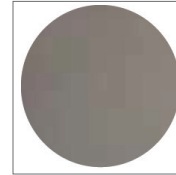


Polyester/Polyester Net Combined

Polyester: Durable, low-friction, easy to clean, apply and remove.

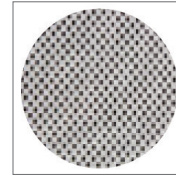


Polyester Net: Soft, thin material that breathes and dries quickly.

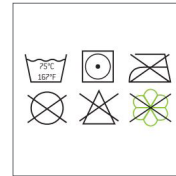

Care

Care instructions located on sling label. Use a hospital grade laundry bag while laundering.

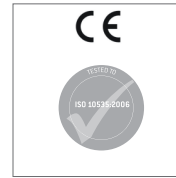

Tested to all relevant standards, including ISO 10535, hoists for the transfer of disabled persons.

Technical Specifications

Part Number	Size	Colour	A	B	C
Polyester/Polyester Net (Combined)					
25120-S	Small		25cm / 9.8"	81cm / 31.9"	84cm / 33.1"
25120-M	Medium		38cm / 15"	98cm / 38.6"	95cm / 37.4"
25120-L	Large		48cm / 18.9"	117cm / 46.1"	105cm / 41.3"
25120-XL	Extra Large		52cm / 20.5"	132cm / 52"	116cm / 45.7"
Maximum User Weight		318kg/701lb			

This product is manufactured in accordance with all relevant ISO standards.



Learn More About This Product



Follow Us on Social Media



facebook.com/humancareworld



instagram.com/humancaregroup



linkedin.com/company/human-care-ab



pinterest.com/humancaregroup



youtube.com/channel/UC7fjInh3zTjRxOP1IDf2Da0

humancaregroup.com