

Bath Sling



Quick Dry
Reinforced Back
Divided Legs



The Human Care Bath Sling can be used in most lifting situations and fits most patients but is most suitable for wet environments. The Polyester Net material quickly dries after use making it ideal for bathing, showering and rehabilitation in pool areas. The Bath Sling designed to give the patient extra support and is a good solution for those users with poor trunk stability.

Quick Dry Material

The Bath Sling is made of a Polyester Net fabric, a soft and thin material that breathes and dries quickly. It is ideal choice for bathing and showering, but can also be used in rehabilitation pools. The sling design creates an upright seating that ensures proper positioning when the patient is placed in a bath chair. Because the material has some elasticity, it can be left behind the patient for a short period of time.

Extra Support

The Bath Sling features a reinforced back and leg parts that offer additional support throughout the back area for those users with poor trunk stability. The Bath Sling can

be used with or without an optional waist belt that is available in an assortment of sizes with a Velcro or buckle closure.

Easy Care and Easy Fit

Our slings are manufactured to the highest quality standards. They feature reinforced straps and high strength nylon webbing around the sling to make them strong enough to resist even the most active use. Each sling label is colour coded to ensure quick and easy selection of the correct sling and the straps are colour coded for easy determination of required positioning.



Accessories



Extension Loops

Used to extend lift straps in order to alter lifting position or to hold leg parts together. Assorted lengths. Part no. 21180.*



Leg Padding and Supports

Ideal for users with sensitive skin or pain. Part no.: 21146-S, M, L Reinforced supports. One size. Part no. 25102



Waist Belts

Velcro or buckle enclosure. Available in Junior and Adult sizing. Velcro Part no. 25145-*** Buckle Part no. 25146-***



Hanger Bars for 2-, 3- or 4-point Suspension

Part no. 57551, 55764 [2- or 3-Point], 55619 [4-Point - Attaches to leg hook]

Lifting Situations



Use to/from a bed

Suitable for transfers and repositioning to and from a bed. Can be left behind patient for short period of time.



Use to/from a chair


Suitable for transfers and repositioning to and from a chair. Can be left behind patient for short period of time.



Use for showering


Only Polyester Net fabric can be used for this application.

Care and Features

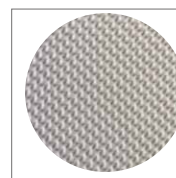


Polyest Net

Soft, thin material that breathes and dries quickly. Ideal for bathing or showering.

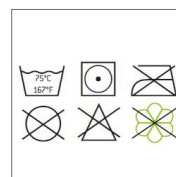

Care


Care instructions located on sling label. Use a hospital grade laundry bag while laundering.

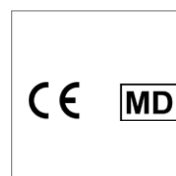

Tested to all relevant standards, including ISO 10535, hoists for the transfer of disabled persons.

Technical Specifications

Part Number	Size	Colour	A	B	C
25035-JS	Junior Small	JS	22cm / 8.7"	56cm / 22"	53cm / 22.9"
25035-JM	Junior Medium	JM	28cm / 11"	63cm / 24.8"	58cm / 22.8"
25035-S	Small	S	35cm / 13.8"	76cm / 29.9"	62cm / 24.4"
25035-M	Medium	M	40cm / 15.7"	86cm / 33.9"	62cm / 24.4"
25035-L	Large	L	44cm / 17.3"	98cm / 38.6"	63cm / 24.8"
25035-XL	Extra Large	XL	48 cm / 18.9"	102cm / 40.2"	64cm / 25.2"
Maximum User Weight	300kg/661lb				

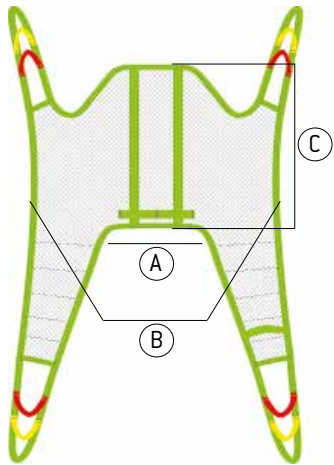

* Available in assorted lengths: 10cm/3.9" (21180-10), 15cm/5.9" (21180-15), 20cm/7.8" (21180-20), 25cm/9.8" (21180-25), 35cm/13.7" (21180-35)

** Junior Sizes: 40cm/15.7" (25145-40B), 50cm/19.7" (25145-50B), 60cm/23.6" (25145-60B); Adult Sizes: 60cm/23.6" (25145-60), 70cm/27.6" (25145-70), 80cm/31.5" (25145-80), 90cm/35.4" (25145-90), 110cm/43.3" (25145-110), 120cm/47.3" (25145-120)

*** Junior Sizes: 40cm/15.7" (25146-40B), 50cm/19.7" (25146-50B), 60cm/23.6" (25146-60B); Adult Sizes: 60cm/23.6" (25146-60), 70cm/27.6" (25146-70), 80cm/31.5" (25146-80), 90cm/35.4" (25146-90), 100cm/39.4" (25146-100), 110cm/43.3" (25146-110), 120cm/47.3" (25146-120)



We improve freedom of movement and quality of care for people

www.humancaregroup.com

HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B
117 43 Stockholm
+46 8 665 35 00
info@humancaregroup.com

HUMAN CARE USA

In8006 Cameron Road, Suite K
Austin, TX 78754 USA
+1 512 476 7199
info.us@humancaregroup.com

HUMAN CARE CANADA

10-155 Colonnade Road
ON K2E 7K1 Ottawa
+1 613 723 6734
info.ca@humancaregroup.com

HUMAN CARE THE NETHERLANDS

Elspeterweg 124
8076 PA Vierhouten
+31 577 41 21 71
info.nl@humancaregroup.com