

# Comfort Sling Soft & Comfort Sling High Soft



Sit on Sling

Undivided Legs

Comfortable



Comfort Sling Soft High  
Art. Nr. 25069



Comfort Sling Soft  
Art. Nr. 25064

Our Comfort Slings are designed to give the patient comfort and a sense of safety. The undivided legs evenly distribute the patient's weight in order to reduce pressure, making this a suitable option for patients with pain. These slings come in two models, one with extra head support.

### Unique Fabric

Our Comfort Slings are made of Soft Air Mesh, a material exclusive to Human Care. This soft material promotes less repositioning making it the optimal choice for patients with pain. In addition, the material does not promote skin breakdown and can therefore be left behind the patient for short periods of time.

### Functional Design

Our Comfort Slings have no reinforcements or handles and feature undivided legs that evenly distribute the patient's weight. With even distribution, pressure is reduced and patient comfort increased. Our Comfort Slings are available in a medium back and high back

model that offers extra head support. Both slings are considered "sit-on" slings and are applied when the patient is lying down.

### Easy Care and Easy Fit

Our slings are manufactured to the highest quality standards. They feature reinforced straps and high strength nylon webbing around the sling to make them strong enough to resist even the most active use. Each sling label is colour coded to ensure quick and easy selection of the correct sling and the straps are colour coded for easy determination of required positioning.



Walking Aids



Healthcare Beds



Lifting Solutions



Hygiene



Convertible Chairs

## Accessories



**Extension Loops**  
Used to extend lift straps in order to alter lifting position or to hold leg parts together. Assorted lengths. Part no.: 21180-\*



**Leg Padding and Supports**  
Ideal for users with sensitive skin or pain. Part no.: 21146-S, M, L Reinforced supports. One size. Part no. 25102



**Upholstered Cover**  
Protects user from hanger bar impact. Not for use in a wet environment. Part no. 50866



**Hanger Bars for 2-, 3- or 4-point Suspension**  
Part no. 57551, 55764 (2- or 3-Point - Quick Release), 55619 (4-Point - Attach to leg hook)


## Lifting Situation




**Use to/from a bed**  
Suitable for transfers and repositioning to and from a bed. Can be left behind patient for short period of time.



**Use to/from a chair**  
Suitable for transfers and repositioning to and from a chair. Can be left behind patient for short period of time.

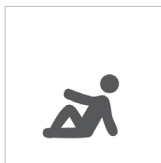


**Use to/from the floor**  
Suitable for transfers from the floor.

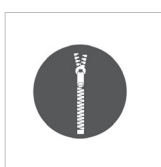


**Use for dressing/undressing**  
Suitable for dressing/undressing applications.

## Care and Features

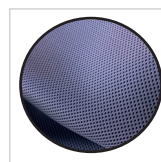

**Soft Air Mesh**  
Lightweight, cool and breathable material. Suitable for sensitive skin.

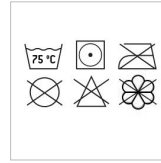

**Care**  
Care instructions located on sling label. Use a hospital grade laundry bag while laundering

## Technical specifications

Part Number	Size	Colour	Model #25064			Model #25069			Diagram	
			A	B	C	A	B	C		
<b>Comfort Sling Soft</b>	<b>Comfort Sing High Soft</b>									
25064-JM	25069-JM	Junior Medium	JM	76cm / 30"	73cm / 28.7"	90cm / 35.4"	76cm / 30"	73cm / 28.7"		111cm / 43.7"
25064-S	25069-S	Small	S	76cm / 30"	80cm / 31.5"	100cm / 39.3"	76cm / 30"	80cm / 31.5"	122cm / 48"	
25064-M	25069-M	Medium	M	86cm / 33.9"	90cm / 35.4"	110cm / 43.3"	86cm / 33.9"	90cm / 35.4"	134cm / 52.8"	
25064-L	25069-L	Large	L	96cm / 37.8"	105cm / 41.3"	116cm / 45.7"	96cm / 37.8"	105cm / 41.3"	139cm / 54.7"	
25064-XL	25069-XL	Extra Large	XL	106cm / 41.7"	118cm / 46.5"	120cm / 47.2"	106cm / 41.7"	118cm / 46.5"	145cm / 57.1"	
Maximum User Weight	300kg / 660 lb									

\* Available in assorted lengths: 10cm/3.9" (21180-10), 15cm/5.9" (21180-15), 20cm/7.8" (21180-20), 25cm/9.8" (21180-25), 35cm/13.7" (21180-35)



*We improve freedom of movement and quality of care for people.*

[www.humancaregroup.com](http://www.humancaregroup.com)



**HUMAN CARE SWEDEN (HQ)**  
Årstaängsvägen 21B  
117 43 Stockholm  
+46 8 665 35 00  
info@humancaregroup.com

**HUMAN CARE GERMANY**  
Industriestraße 13, Gebäude 1/9  
63755 Alzenau  
+49 (0) 6023 994 936 8  
infode@humancaregroup.com

**HUMAN CARE CANADA**  
10-155 Colonnade Road  
ON K2E 7K1 Ottawa  
+1 613 723 6734  
info.ca@humancaregroup.com

**HUMAN CARE THE NETHERLANDS**  
Elspeterweg 124  
8076 PA Vierhouten  
+31 577 41 21 71  
info.nl@humancaregroup.com

**HUMAN CARE UNITED STATES**  
8006 Cameron rd, Suite K  
TX 78754 Austin  
+1 5 124 76 7199  
info.us@humancaregroup.com