



Easy to Store
Homecare/LTC
Wood Option



Our Floorline LTC is the optimal bed solution for home care use or long-term care facilities. Its one piece design makes it easy to fold, store and transport. The Floorline LTC comes in two convenient options including one with a wood-grain finish that compliments any room.

Clinical and Patient Benefits

The Floorline LTC lowers to a distance of 9.5cm/3.7" from the floor. This distance, although seemingly small, has enormous benefits to both the patient and the facility. By just adding 5cm to the distance between the bed and the floor, the impact force of a fall is increased by 50%. This can be the difference between a simple bruise or a fractured hip for the patient that may result in an extended hospital stay or even death. In addition to the physical effects, fall injuries also have psychological effects on the patient as they experience a loss of confidence and fear after a fall.

Optimal Features

The Floorline LTC features standard positions

normally found only on hospital beds, such as Trendelenburg, reverse Trendelenburg and cardiac positioning. With a full assortment of accessories, an easy-to-use handset, and durable LINAK electronics, the Floorline LTC is the best choice for patient care.

Quality

Our products are designed and manufactured to meet all the requirements of relevant standards. At every stage of the development process, every single product is given extensive attention in design, development and manufacturing. This means that a user of this product feels confident with its reliability and the caregiver feels confident with its safety.



Walking Aids



Healthcare Beds



Lifting Solutions




Hygiene



Convertible chairs

Accessories

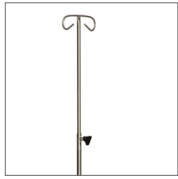

IV Pole
Adjustable height. Can be mounted on either end of bed. Maximum SWL 7kg/15.4lb. Part no. IVP



Care Assist Rail
Designed to assist patients safely in and out of a bed. Ergonomic design. Part no. 90702



Extension Kit
Extends bed by 17.5cm/7" for taller patients. Part no. 90709



Handset
Five function - included with bed. Part no. 75746
Handset holder Part no. LC2250, Lock Out Key Part no. 75790



Self Help Pole
Assists patient moving within confines of bed. Maximum SWL 75kg/165.3lb. Part no. SHP600



Side Rails
1/4 size, art. nr: 72794
1/2 size, art. nr: 72795
3/4 size, art. nr: 72797



Wall Bumper Bar
Protects the head of the bed. Mounted on lower support bracket between wheels. Part no. 90705



Oxygen Holder
Carries C size tank. Attaches to legs. Part no. 90707

Features



Compact
Folds compactly for storage or transport.



Standard Positions
Includes auto-contour, backrest, kneebreak, Trendelenburg, reverse Trendelenburg, cardiac.



Premium Casters
One central locking and three steering/directional casters for premium maneuverability.



Manufactured to all relevant standards including IEC 60601-2-52:2009, EN1970:2002, AS/NZS 3200.2.38:2007

Technical Specifications

Floorline LTC			
Part no. White Finish	90600 (EU), 90602 (UK), 90606 (NA), 90608 (AU)	Part no. Woodgrain Finish	90601 (EU), 90603 (UK), 90607 (NA), 90609 (AU)
Platform Position (Min-Max)	9.5cm - 70cm / 3.7" - 27.5"	Overall Bed Length (Standard to Extended with Extension Kit)	215cm - 232.5cm / 84.5" - 91.5"
Platform Dimensions (L x W)	215cm x 90.5cm / 84.5" x 35.5"	Platform Dimensions with Extension Kit (L x W)	232.5cm x 90.5cm / 91.5" x 35.5"
Sleeping Dimensions (L x W)	200cm x 90cm / 78.7" x 35.5"	Sleeping Dimensions with Extension Kit (L x W)	217.5cm x 90cm / 85.5" x 35.5"
Output Voltage	Max. 24V DC	IP Rating	IPX4
Input Voltage	120V, 60Hz (90404 to 90607), else 230V, 50Hz	Product Weight	106kg / 234lb
Platform Panel Angles	Backrest 70°, Thigh 45°, Knee 110°, Calf 25°, Trendelenburg/Reverse Trendelenburg 18°	Safe Working Load (SWL)	200kg / 440lb
Duty Cycle	2 minutes operation, 18 minutes rest	Maximum User Weight	165kg / 363lb

This product is manufactured in accordance with all relevant ISO standards.

