



Stable
Compact design
Ergonomic



Electric, height adjustable tilting commode/shower chair

Our Dolphin II allows personal hygiene activities to be carried out in the most ergonomic manner. While the client is sitting in a relaxed, tilted position, the optimal height adjustment function prevents the caregiver from having to stoop.



The large range of the seat height means that the caregiver can even address the personal hygiene needs of the client at hip height. The seat and backrest tilt as one, keeping the client in a stable position. The electrically operated tilt is seamless. It is easy for the caregiver to operate the system using only one hand and without any effort, thus facilitating the most comfortable position for both the client and caregiver. It is easier for the client to stand up from the chair from a sitting position because of the forward tilt.

Dolphin II has a headrest that can be adjusted in terms of height, depth and angle. The bowl-shaped design of the headrest means that the client can rest their head securely and comfortably on the headrest.

The level of the footplate can be adjusted to provide correct foot support for clients of different heights.



The armrests continue around the front so that the client has a secure support, which is useful when it is necessary to lean forward to receive personal care. This design also makes the client feel safe. The armrests are there when needed, but can be easily folded up so that they do not get in the way of a transfer in/out of the chair.

Due to its compact size the Dolphin II is ideal for use in small spaces, both at home, narrow doorways and toilets/bathrooms as well as in healthcare institutions where other necessary medical devices can often clutter an area. The chair is easy to maneuver thanks to the large, high-quality twin castors. This, and many other factors, make Dolphin II, a durable quality product.



Accessories



Toilet bucket
Toilet bucket, including lid
Article number
04-M0312-0016



Toilet seat
Toilet seat with lid for the keyhole shaped opening for extra comfort.
Article number
04-M0307-0009



Battery
Extra Linak battery to ensure the lift is always ready for use.
Article number
04-N0101-0012



External charger
Charger for external battery charging.
Article number
04-N0105-0018


Use



Personal hygiene
Provides support to clients for their personal care, such as intimate hygiene and getting dressed and undressed.



Ergonomic
Excellent maneuverability and height adjustment will improve the caregiver's posture, reducing the strain on their body.

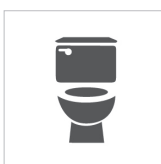

Toileting
Can be wheeled into toilet area. Provides support during toileting. Can be positioned over a suspended toilet. Can be used with a toilet bucket.



Showering
Can be wheeled into the bathroom. Offers comfortable and safe support during showering and washing (IPX5, waterproof).


Features



Casters
4 twin castors, with brakes (Ø 125 mm).



Hand control
User-friendly hand control, easy to use, ergonomic design with durable coiled cable.

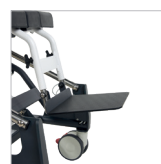

Footplate
Height adjustable, horizontal anti-slip footplate to keep feet in place while the client is toileting. Can be folded away during transfers.



Headrest
Height, depth and angle-adjustable, bowl-shaped PUR headrest for safe and comfortable support of the neck and head.

Technical specifications

Specifications	Dolphin II, 90561
Height	920 mm / 36.2"
Width	670 mm / 26.3"
Seat height	480-1000 mm / 18.8-39.3"
Position of the push handle (min/max)	1020-1310 mm / 40.1-51.5"
Size of the seat	420x430 mm / 16.5x 16.9"
Tilt	-5° to +26°
Diameter of castors with brake	125 mm / 4.9"
Maximum user weight	180 kg / 396lbs
Maximum weight per arm rest	80 kg / 176lb
Turning radius	1100 mm / 43.3"
Water resistance	IPX5 (waterproof)
Water resistance of hand control and lift actuator	IPX6

Specifications	Dolphin II, 90561
Headrest	height, depth and angle-adjustable
Rechargeable battery (incl. charging cable)	Capacity 2.9 Ah
Charger, integrated in control box	100-240/50/60 V/HZ
Material of the frame	Powder-coated steel
Material of seat, backrest and headrest	PUR (polyurethane)
Toilet seat	with lid and integrated guide rails for toilet bucket
Armrests	individually foldable, with integrated front support
Footplate	foldable, height adjustable, with anti-slip coating
Hand control functions	Up/down, tilting, battery status, service indicator
Accessory	Safety belt (90562), adjustable length

99494_Rev 01 ENG Product Sheet Dolphin II, Feb 2021



We improve freedom of movement and quality of care for people.

www.humancaregroup.com



HUMAN CARE SWEDEN (HQ)
Årstaängsvägen 21B
117 43 Stockholm
+46 8 665 35 00
info@humancaregroup.com

HUMAN CARE GERMANY
Industriestraße 13, Gebäude 1/9
63755 Alzenau
+49 (0) 6023 994 936 8
infode@humancaregroup.com

HUMAN CARE CANADA
10-155 Colonnade Road
ON K2E 7K1 Ottawa
+1 613 723 6734
info.ca@humancaregroup.com

HUMAN CARE THE NETHERLANDS
Elspeterweg 124
8076 PA Vierhouten
+31 577 41 21 71
info.nl@humancaregroup.com

HUMAN CARE UNITED STATES
8006 Cameron rd, Suite K
TX 78754 Austin
+1 5 124 76 7199
info.us@humancaregroup.com