




Fitting Protocol

Back Belt

Models 25130, NSB-500/501/501C

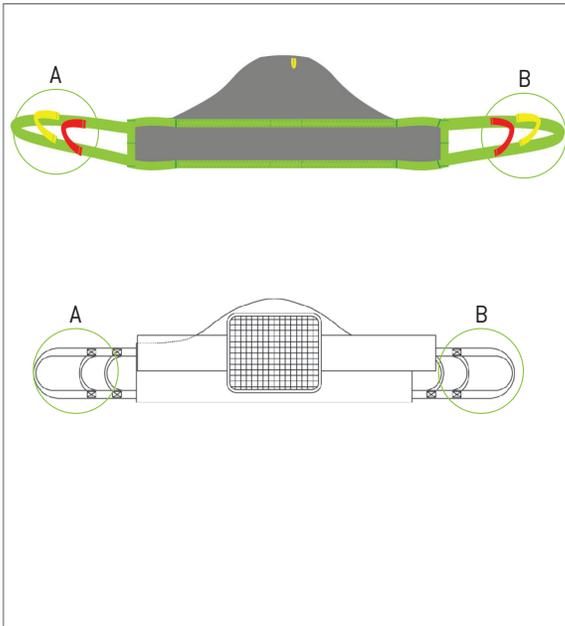

Strap A

- Red Strap
- Yellow Strap
- Green Strap

Strap B

- Red Strap
- Yellow Strap
- Green Strap

- Model 25130 - Polyester**
- Model NSB-500 - Polyester**
- Model NSB-501 - Cotton, Single Patient Use**
- Model NSB-501C - Cotton, Single Patient Use, Clips**

Choose the loop combination that will give the best distance from the Sit-to-Stand for the intended use. When securing straps to Sit-to-Stand device, match colours on both sides of the patient.

Maximum User Weight: 150kg/330.7lb

Transfer Situation:

- Raise/Lower To/From a Bed
- Raise/Lower To/From a Chair
- Gait Training/Mobilization
- Raise/Lower To/From a Toilet

Available sizes:

- Indicate size required
- S
 - M
 - L
 - XL
- Note: NSB-501 size XL is not available

Leg Placement:

- Not applicable for this sling
Divided relaxed position
- Not applicable for this sling
Open position
- Not applicable for this sling
Divided position
- Not applicable for this sling
Closed position

Suspension:

- 2-Point

Accessories required:

Patient Information:

Name	
Facility	Room
Lift Used	Hanger Bar Used

HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B
117 43 Stockholm
Phone:- +46 8 665 35 00
Fax:- +46 8 665 35 10
info.se@humancaregroup.com

HUMAN CARE AUSTRALIA

Unit 4, 49 Bertie Street
Port Melbourne, VIC 3207
Phone : +61 1300 132 736
Fax : + 61 (03) 9645 4954
info.au@humancaregroup.com

HUMAN CARE CANADA

10-155 Colonnade Road
Ottawa, ON K2E 7K1
Phone: +613.723.6734
Fax: +613.723.1058
info.ca@humancaregroup.com

HUMAN CARE NEDERLAND

Elspeterweg 124
8076 PA Vierhouten
Phone: +31 577 412 171
Fax: +31 577 412 170
info.nl@humancaregroup.com

HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K
Austin, TX 78754
Phone: +512.476.7199
Fax: +512.476.7190
info.us@humancaregroup.com

Back belts are designed for use with a Sit-to-Stand device and help safely raise or lower patients from a bed, chair or toilet. The wide reinforced back support allows for gait training or mobilization exercises when the float plate is removed from the device. Back belts are available in reusable Polyester models or single patient use cotton models. These belts are recommended for patients with poor trunk stability and difficulty standing. Some leg control is necessary to use this product.

Additional Instructions

To staff/caregiver:

To patient:

Fitting Protocol Completed by:

Name	
Date	Position/Title



Important! This fitting protocol does not replace the user manual. Read all instructions before use!

Additional Documentation

Can be found on our website at www.humancaregroup.com/product-category/slings



Quick Reference Guide



User Manual



Sling Care Instructions



Sling Inspection Form

