

General Purpose Sling



Single Patient Use
NEW Steering handle
Loop or Clip



Our General Purpose slings have a high weight capacity and a generous overall width to offer a stretched seating position that is well-suited for lifting heavier users in all situations. The integrated head support allows the patient to be lifted even from a floor.



Promotes Infection Control

These single patient use slings offer the convenient option of disposal after use combined with the added benefit of hygiene control. Since the sling can only be used by one patient at a time, facility infection control protocols are easily supported.

Two Convenient Options

The General Purpose sling is available with either reinforced loops that ensure durability and ease of use, or with clips for easy and secure attachment. Either option is an ideal solution for sudden positioning change requirements. These slings are available in adult sizes Small to Extra Large and can carry a maximum capacity of 318kg/700lb.

Optimal Features

These slings feature divided legs and a high back for added support. They are made of a strong but soft polyester/cotton blend and are ideal for a seated, horizontal or semi-recumbent transfer. The unique material wicks away fluids and can be left behind a patient. Model #NSB-104 can be used with a spreader bar. Model #NSB-105 features loops that can be attached to a hanger bar with 2-point, 3-point or 4-point suspension.



Accessories



Extension Loops

Used to extend lift straps in order to alter lifting position or to hold leg parts together. Assorted lengths. Part no.: 21180-*



Leg Padding

Soft padded supports for legs. Ideal for users with sensitive skin or pain. Synthetic sheepskin. Part no.: 21146-S, M, L



Hanger Bar

For NSB-105. For 2-, 3-, or 4-point suspension. Part no. 57551, 57664 [2- or 3-point], 55619 [4-point - attach to leg hook on lift]



Spreader Bar

For NSB-104. Powered Spreader Bar Part no.: Small 92441, Medium 92442, Large 92443

Lifting Situations



Use to/from a bed

Suitable for transfers and repositioning to and from a bed. Can be left behind patient for short period of time.



Use to/from a chair

Suitable for transfers and repositioning to and from a chair. Can be left behind patient for short period of time.



Use to/from the floor

Suitable for transfers from the floor.



Single Patient Use

Sling is intended for use by a single patient at a time. Dispose after use.

Care and Features




Clip or Loop Sling

Two convenient attachment options - clips or loop slings.



Reinforced Straps

Colour coded reinforced straps for easy determination of required positioning.

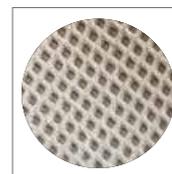


Non-woven polypropylene

Material allows fluids to wick away from patient. Do not launder. Single Use.



Technical Specifications

Article number	Size	Color	A	B	C
NSB-104C-S	Small	S	30cm / 12"	84cm / 33"	78cm / 31"
NSB-104C-M	Medium	M	43cm / 17"	88cm / 35"	97cm / 38"
NSB-104C-L	Large	L	53cm / 21"	103cm / 41"	114cm / 45"
NSB-104C-XL	Extra Large	XL	68cm / 27"	114cm / 45"	133cm / 53"
NSB-105-S NSB-105V-S ¹	Small	S	30cm / 12"	84cm / 33"	78cm / 31"
NSB-105-M NSB-105V-M ¹	Medium	M	43cm / 17"	88cm / 35"	97cm / 38"
NSB-105-L NSB-105V-L ¹	Large	L	53cm / 21"	103cm / 41"	114cm / 45"
NSB-105-XL NSB-105V-XL ¹	Extra Large	XL	68cm / 27"	114cm / 45"	133cm / 53"
Maximum user weight	318kg/700lb.				

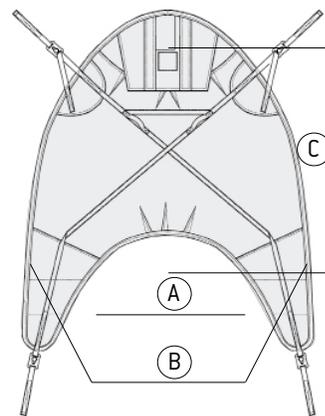

* Available in more sizes: 10cm/3.9" [21180-10], 15cm/5.9" [21180-15], 20cm/7.8" [21180-20], 25cm/9.8" [21180-25], 35cm/13.7" [21180-35]

www.humancaregroup.com

HUMAN CARE SWEDEN (HQ)
Årstaängsvägen 21B
117 43 Stockholm
+46 8 665 35 00
info@humancaregroup.com

HUMAN CARE CANADA
10-155 Colonnade Road
ON K2E 7K1 Ottawa
+1 613 723 6734
info.ca@humancaregroup.com

HUMAN CARE THE NETHERLANDS
Elspeterweg 124
8076 PA Vierhouten
+31 577 41 21 71
info.nl@humancaregroup.com

HUMAN CARE UNITED STATES
8006 Cameron rd, Suite K
TX 78754 Austin
+1 5 124 76 7199
info.us@humancaregroup.com