

# General Purpose Sling



Reusable  
Padded Support  
Washable



Our General Purpose sling is a reusable sling that has a high weight capacity and a generous overall width to offer a stretched seating position that is well-suited for lifting heavier users in all situations. The integrated head support allows the patient to be lifted even from a floor and its unique material can be laundered to allow for use by multiple patients.

#### Ventilated Polyester Material

The NSB General Purpose Sling is a versatile sling that can be used for many different patients with many different lifting situations. It has a padded head support and is made from a polyester material with built-in flexibility to achieve proper leg and body positioning. This flex material allows the sling to remain behind the patient for short periods of time and the ventilated portion allows breathability.

#### Bariatric Option

The NSB General Purpose sling is available in adult sizes Small to Extra Large and can carry a maximum capacity of 318kg/700lb as well as a bariatric version that can carry

a maximum capacity of 454kg/1000lb. This product is only available in North America.

#### Optimal Features

This sling features divided and padded legs and a high back with an integrated padded head support for added comfort. The strong, durable material and design makes it ideal for a seated or semi-recumbent transfer. The sling features 8 straps for optimal positioning on both a ceiling or floor lift. The unique material allows the sling to be laundered for use by multiple patients.

**Note: This product is only sold in North America.**

#### HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B  
117 43 Stockholm  
Phone:- +46 8 665 35 00  
Fax:- +46 8 665 35 10  
info.se@humancaregroup.com

#### HUMAN CARE CANADA

10-155 Colonnade Road  
Ottawa, ON K2E 7K1  
Phone: +613.723.6734  
Fax: +613.723.1058  
info.ca@humancaregroup.com

#### HUMAN CARE NEDERLAND

Elspeterweg 124  
8076 PA Vierhouten  
Phone: +31 577 412 171  
Fax: +31 577 412 170  
info.nl@humancaregroup.com

#### HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K  
Austin, TX 78754  
Phone: +512.476.7199  
Fax: +512.476.7190  
info.us@humancaregroup.com

## Accessories



### Extension Loops

Used to extend lift straps in order to alter lifting position or to hold leg parts together. Assorted lengths. Part no.: 21180-\*



### Leg Padding

Soft padded supports for legs. Ideal for users with sensitive skin or pain. Synthetic sheepskin. Part no.: 21146-S, M, L



### Hanger Bar

For 2-, 3-, or 4-point suspension. Part no. 57551, 57664 (2- or 3-point), 55619 (4-point - attach to leg hook on lift)



### Ceiling or Floor Lifts

Can be used with most standard floor lifts. Compatible with Altair, Roomer S or HeliQ lifts.

## Lifting Situations



### Use to/from a bed

Suitable for transfers and repositioning to and from a bed. Can be left behind patient for short period of time.



### Use to/from a chair

Suitable for transfers and repositioning to and from a chair. Can be left behind patient for short period of time.



### Use to/from the floor

Suitable for transfers from the floor.



Tested to all relevant standards, including ISO 10535, hoists for the transfer of disabled persons.

## Care and Features




### Padded Legs and Head

Extra padding for optimal comfort and patient positioning.




### Reinforced Straps

Colour coded reinforced straps for easy determination of required positioning.

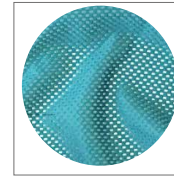

### Ventilated Polyester

Material allows fluids to wick away from patient.

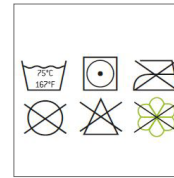


### Care


Care instructions are located on the sling label. Use a hospital grade laundry bag while laundering.

## Technical Specifications

Part Number	Size	Colour	A	B	C
NSB-103-S	Small	S	23.5cm / 9.25"	77.5cm / 30.5"	78cm / 31"
NSB-103-M	Medium	M	42.5cm / 16.75"	97cm / 38.25"	97cm / 38"
NSB-103-L	Large	L	53cm / 21"	114cm / 45"	114cm / 45"
NSB-103-XL	Extra Large	XL	68cm / 27"	133cm / 52.5"	133cm / 53"
NSB-109-2XPH	Bariatric	XXXL	62.2cm / 24.5"	144.5cm / 57"	85.1cm / 33.5"
NSB-103C-S	Small	S	30cm / 12"	84cm / 33"	78cm / 31"
NSB-103C-M	Medium	M	43cm / 17"	88cm / 35"	97cm / 38"
NSB-103C-L	Large	L	53cm / 21"	103cm / 41"	114cm / 45"
NSB-103C-XL	Extra Large	XL	68cm / 27"	114cm / 45"	133cm / 53"
Maximum User Weight	Small to X Large: 318kg / 700 lb Bariatric: 454kg / 1000lb				

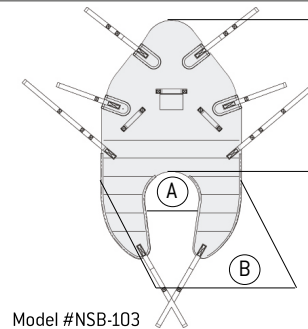

Model #NSB-103

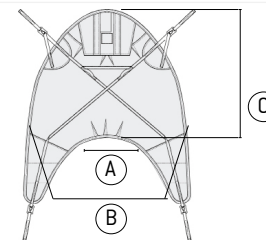

Model #NSB-103C

\* Available in assorted lengths:  
10cm/3.9" (21180-10),  
15cm/5.9" (21180-15),  
20cm/7.8" (21180-20),  
25cm/9.8" (21180-25),  
35cm/13.7" (21180-35)