Levanto











The Levanto is a safe and comfortable standing and transfer device designed to help patients with some trunk balance and weakened leg muscles to stand up with a natural motion pattern. It features clip sling attachments and its compact design makes it the ideal option for smaller spaces.

Optimal Features

The Levanto is an electric sit-to-stand device that encourages a natural standing motion and is designed for standing, sitting, dressing, toileting and balance training for users with low core or leg strength. The shin pad angle adjusts with the patient while they are being raised for optimal comfort and security and the legs are easily spread electronically to allow direct access while seated (for example, in a wheelchair). The non-slip surface footplate offers stable support and is removable for easier cleaning and hygiene control.

Ideal for Home Care Use

The Levanto features a compact footprint that allows it to easily pass through narrow

doorways. It features a lifting arm with additional handles for comfort and easy usage. The shorter legs allow for easier rotation, especially in a home environment that may have less space for manoeuvring. With a LINAK battery and actuator, you are ensured years of trouble-free use.

Quality

Our products are designed and manufactured to meet all the requirements of relevant ISO, UL, IEC and EN standards. At every stage of the development process, every single product is given extensive attention in design, development and manufacturing. This means that a user of this product feels confident with its reliability and the caregiver feels confident with its safety.

HUMAN CARE SWEDEN (HQ) Årstaängsvägen 21B 117 43 Stockholm Phone:- +46 8 665 35 00 Fax:- +46 8 665 35 10 info.se@humancaregroup.com

HUMAN CARE AUSTRALIA Unit 4, 49 Bertie Street Port Melbourne, VIC 3207 Phone : +61 1300 132 736 Fax : + 61 (03) 9645 4954 info.au@humancaregroup.com

HUMAN CARE CANADA

10-155 Colonnade Road Ottawa, ON K2E 7K1 Phone: +613.723.6734 Fax: +613.723.1058 info.ca@humancaregroup.com HUMAN CARE NEDERLAND Elspeterweg 124 8076 PA Vierhouten Phone: +31 577 412 171 Fax: +31 577 412 170

info.nl@humancaregroup.com

HUMAN CARE UNITED STATES 8006 Cameron Road, Suite K Austin, TX 78754 Phone: +512.476.7199 Fax: +512.476.7190 info.us@humancaregroup.com

Accessories



Back Belts Assortment of sizes and designs to accommodate most every lifting need.



Battery Complete system with control box enclosed in a discreet module. Included with lift.



Hand Control Slow or fast setting for arm up/down, electrical leg spreading. Included with lift.



Levanto

Wheels Rear wheels feature kick-

Technical Specifications


brake. Diameter 125mm.

Lifting Situations



Use to/from a bed Use to raise/lower a patient from/to a bed.

Use to/from a chair Use to raise/lower a patient from/to a chair.

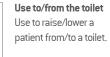



Control Box

Emergency stop, emergency lowering and soft-start functions.



Lifting Arm Lift arm extends 80cm / 31.5" to accommodate taller patients.




Spreadable Base Electrical adjustment.



Tested to all relevent standards, including ISO 10535, hoists for the transfer of disabled persons.

Loranto				
Part No.	20000170	IP Rating	IPX4	
Width	68.7cm / 27"	Turning Radius	128cm / 50.3"	
ength	110cm / 43.3"	Front Wheels Rear Wheels with Brake	10cm/4" 12.5cm/5"	
leight	Arm Lowered: 88.9cm / 35" Arm Raised ¹ : 164.8cm / 64.9"	Battery	24V DC/2.9Ah, 2x 12V/2.9Ah	11137 -
leight from Floor	11.9cm/4.7"	Weight	49kg/108lb	530
djustable Leg Min - Max)	Outside: 67.9cm - 113.1cm /26.7" - 44.5"	Weight Capacity	160kg/352.7lb	
Sound Level	48dB (loaded) while lifting 52dB (loaded) while spreading legs			

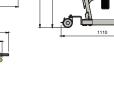

¹Measurement based on arm mast at highest setting

This product is manufactured in accordance with all relevant ISO standards



Learn More About This Product





Follow Us on Social Media

facebook.com/humancarecda

0 instagram.com/humancaregroup





pinterest.com/humancaregroup

youtube.com/channel/UC7fyiNh3zTjRx0P1IDf2DaQ

You Tube



humancaregroup.com